



Spider vein removal machine

LB-100



PRODUCTS DISPLAY



Packaging

PRODUCT ADVANTAGES

✔ 30MHz high frequency, painless & effective

✔ See results in only one treatment

✔ Not affect normal work and life

✔ Security by Micro Control Technology


✔ Touchscreen and self-checking system

✔ Gold quality with 3 years warranty


✔ Fast Reinvest, minimum consumption cost




HANDLE DETAILS

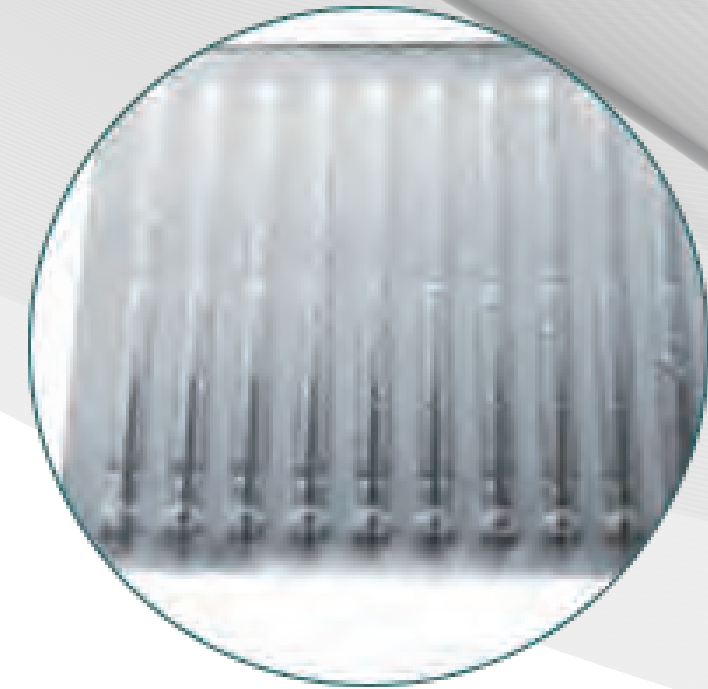

NEEDLE DETAILS

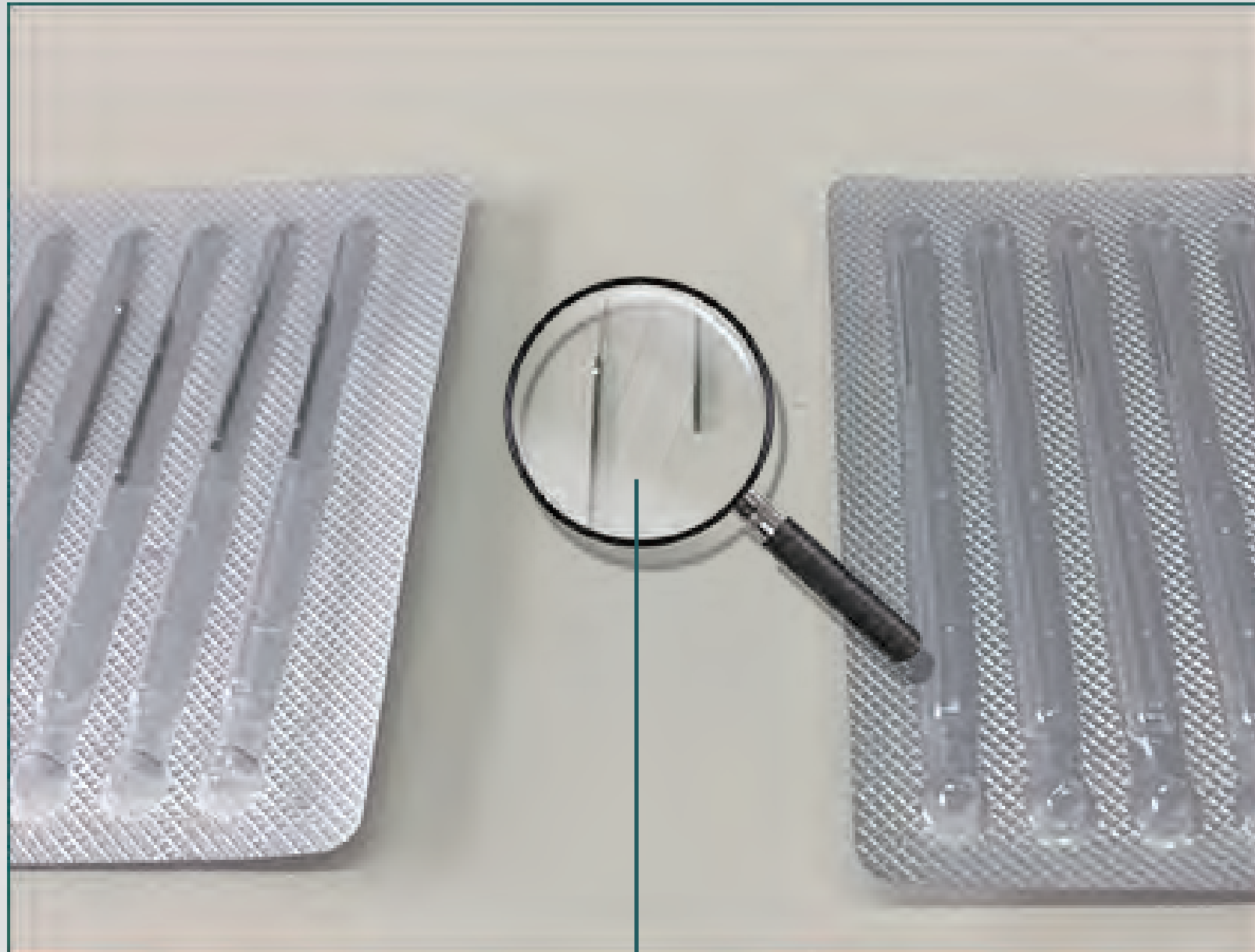

0.03mm needle for skin tag-
es, pigmentation removal

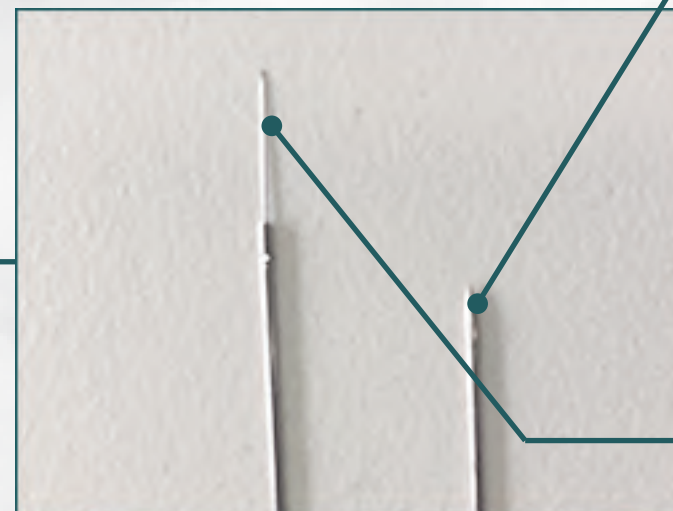



0.01mm needle for spider
veins removal

PRODUCT PARAMETERS


Model	LB-100
Output frequency	30MHz
Output power	0-100% (adjustable) $\pm 10\%$
Input power	150W
Output voltage	AC 30-150V
Treatment needle	0.01mm(diameter)
Output mode	Continuous and pulsed two modes
Liquid Crystal Display	8.4 Inch
Input voltage	220V-50HZ/ 110V-60HZ
Operating mode	Touch

WORKING PRINCIPLE


This equipment applies painless high frequency technology, and can make a thermal coagulation of Hemoglobin and endothelial cells inside the dilated capillaries in a millisecond to close the dilated capillaries immediately. The diameter of its treatment needle is only 0.01mm, so the damage to skin epidermis could be omitted when receiving treatments. The extremely high energy density of this apparatus can do a coagulation to the dilated capillaries instantly, and nearly any skin tissues or coagulated tissues will fall off in 3-7 days, so this equipment can also effectively treat pigmented lesions, like proliferative lesions, age spots, sunburn etc.

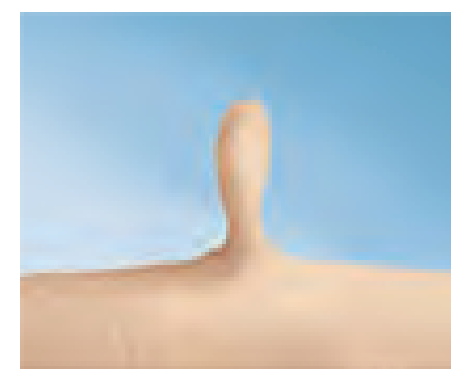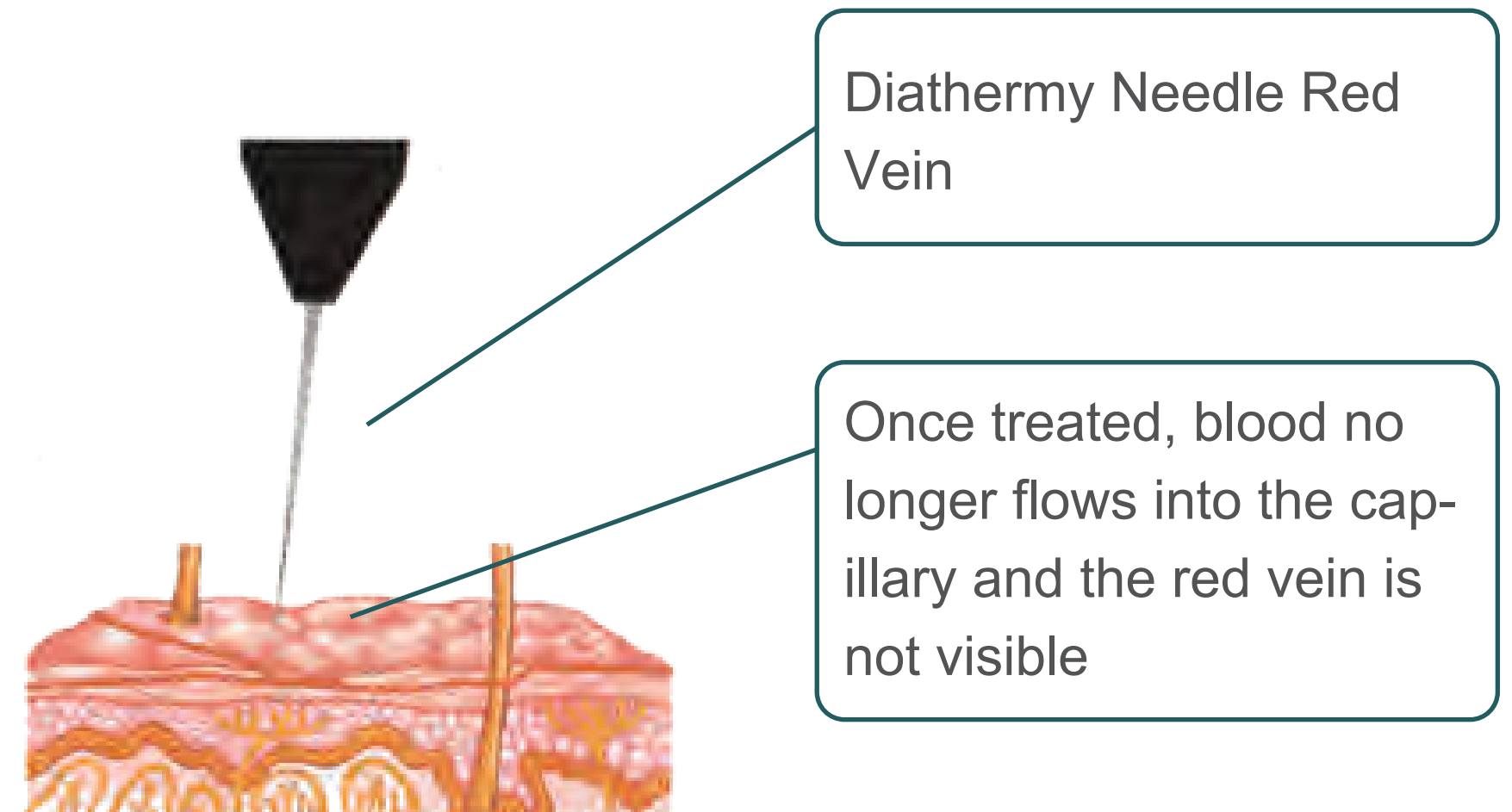

High frequency energy acts instantly to thermal coagulate the skin tag from the base.

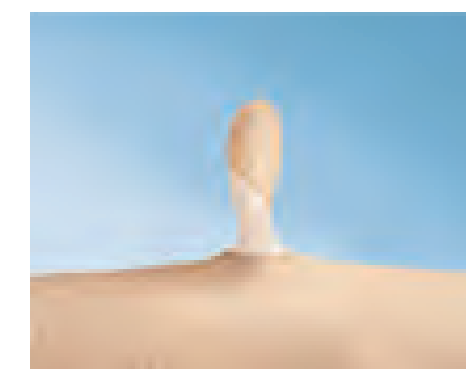


Skin tag base turns white within second as tissue is being destroyed

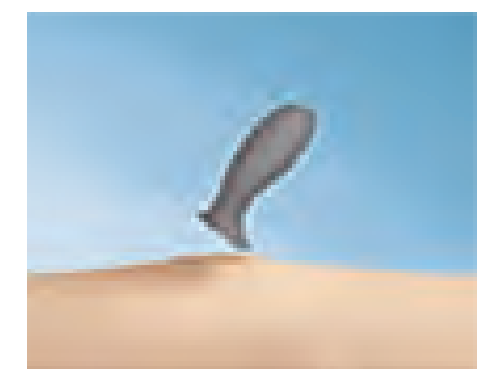





Skin tag falls off after new heal the skin has formed under the treated area in 3-7 days

TREATMENT SCOPE

1. All kinds of redness, such as spider-web like veins, linear redness, filamentous redness, genetic redness, naturally occur redness, flat cherry-like vascular tumors.
2. Benign lesions, such as skin bumps, milia, flat wart, handover nevus, compound nevus, fat particles and so on; also pigmentary lesions, such as, spots, age spots, sun spots, pigmentation and so on.

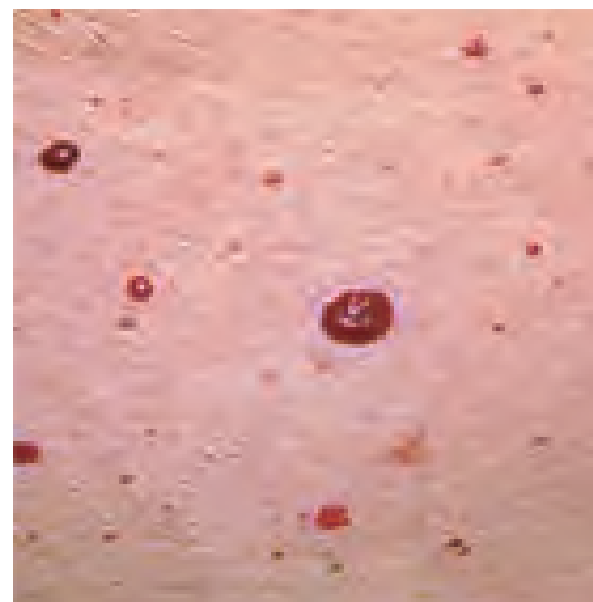





BEFORE VS AFTER

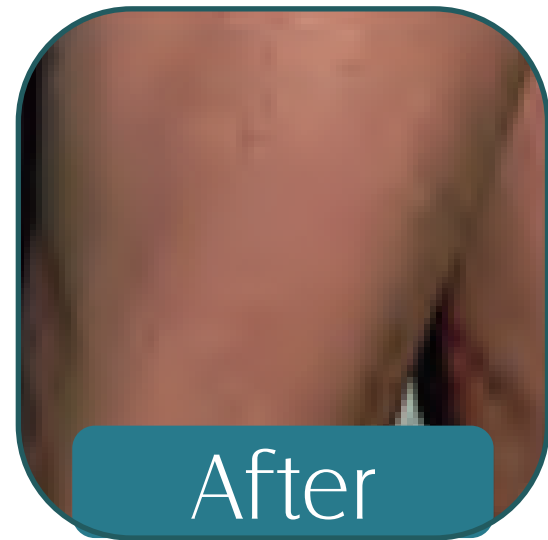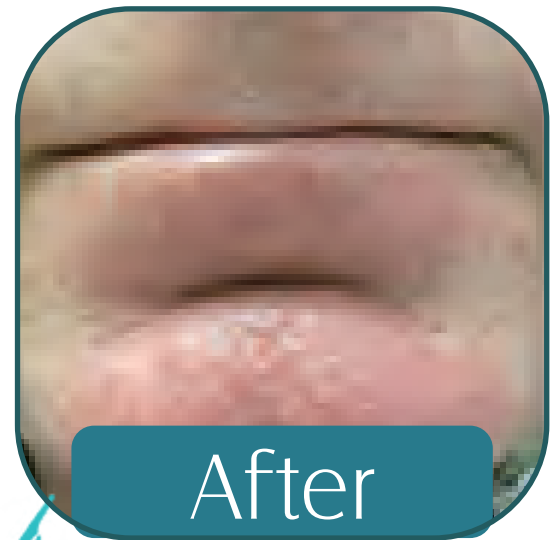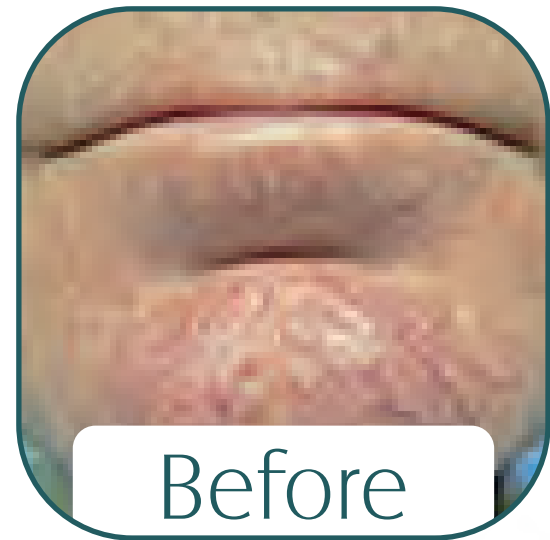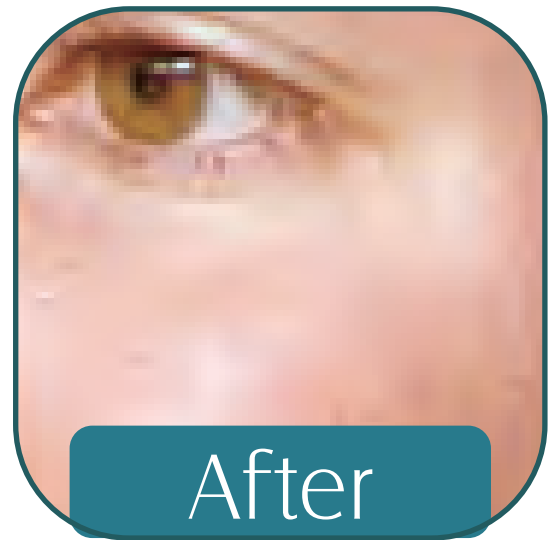

VS



CUSTOMER FEEDBACK

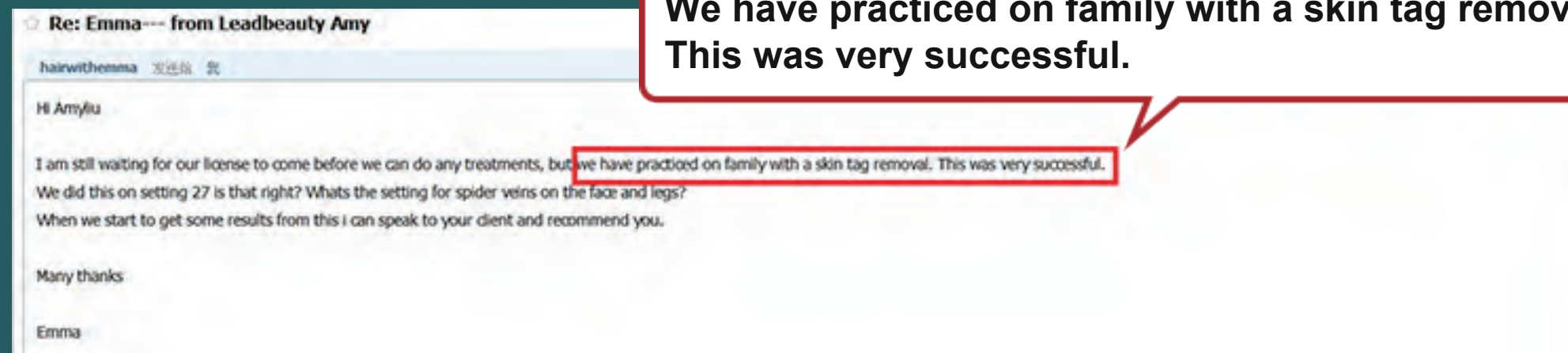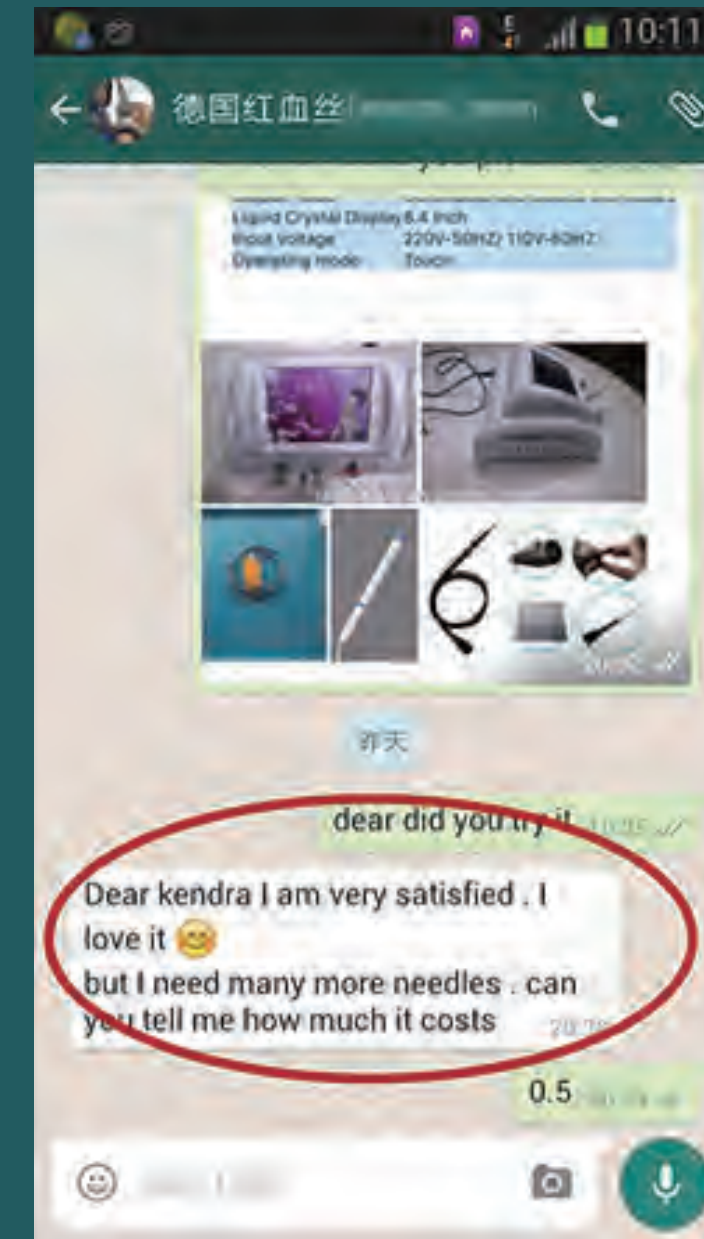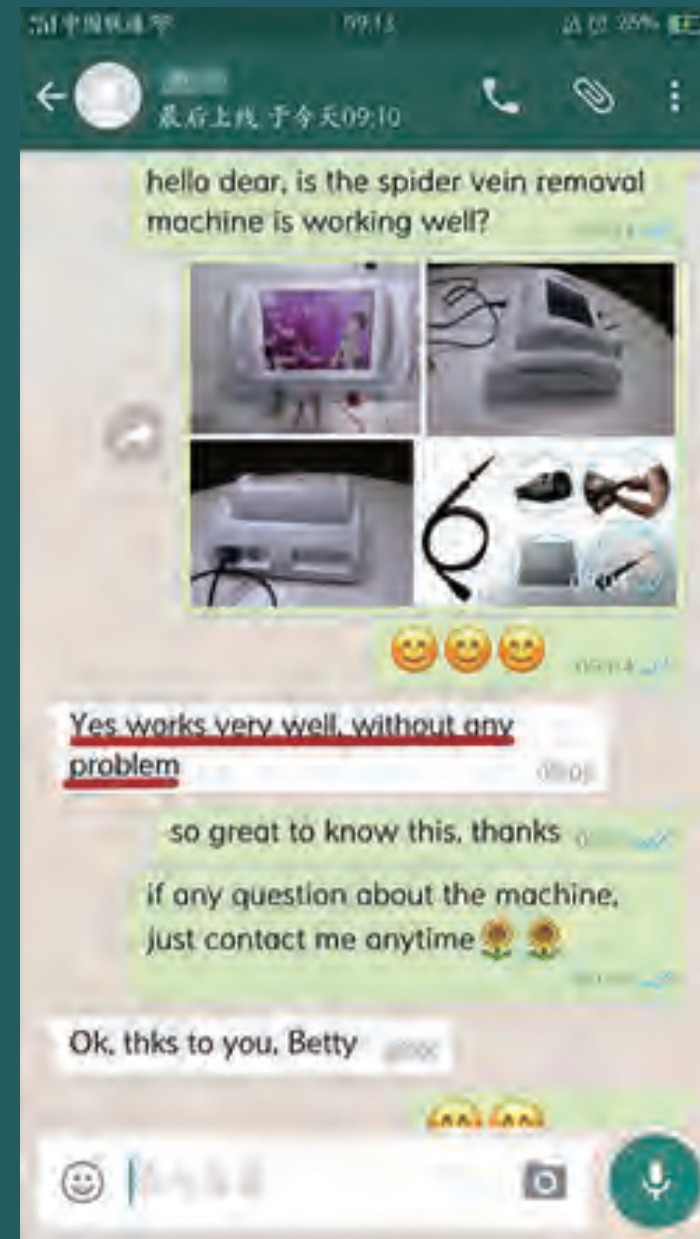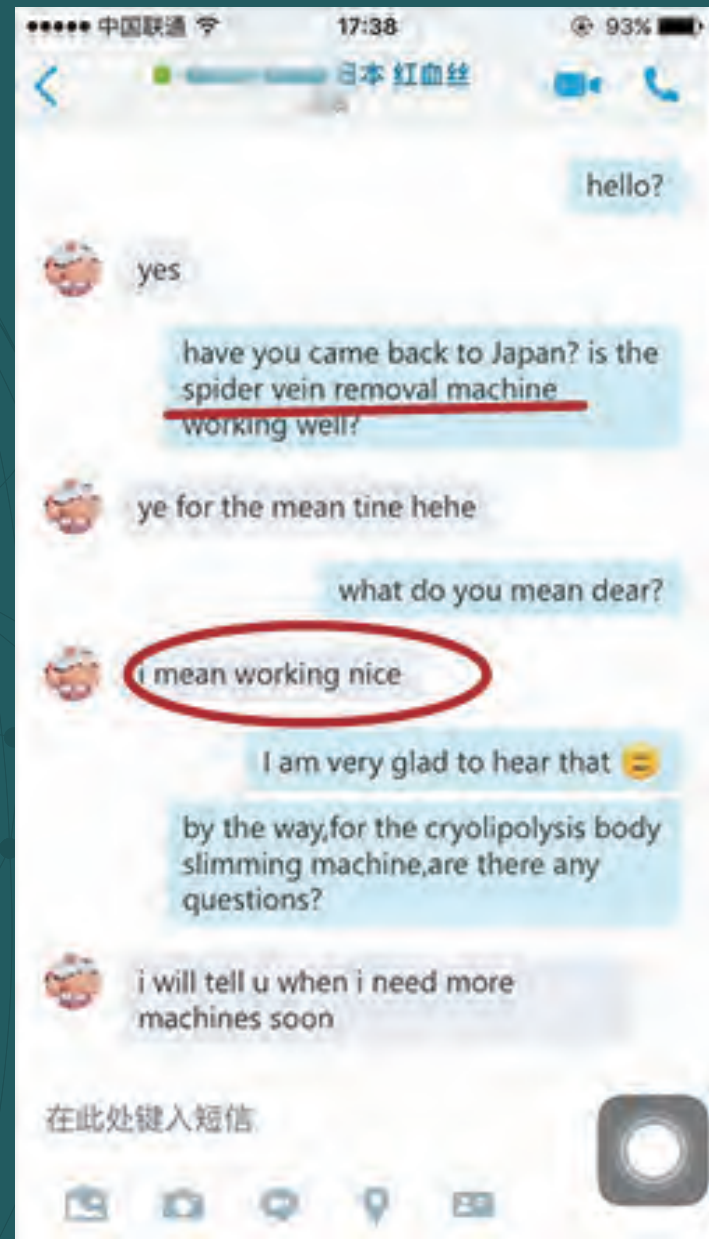

We have practiced on family with a skin tag removal. This was very successful.

About us



Over 13 years manufacturer- Strong production capacity Beijing Leadbeauty International S&T CO. Ltd can keep you always ahead in the aesthetics & medical devices field. LEADBEAUTY has its own research & development, sales and after-sales departments; can offer the professional technology supports at the first time. We are one of the most leading & reliable manufacturers in China, have a professional team integrating with optics, machinery and electricity that keep us ahead on in this field.



Beijing LEADBEAUTY International S&T CO. Ltd

Tel: 86-10-60418895

Mob & whatsapp: +86 13269056123

E-mail: leadbeauty@leadbeauty.net

Website: www.leadbeauty.net

Addr: West zone No.5 Mauhwa Industry Park, No.1 Caida 3rd Street, Shunyi District, Beijing, China